

The "Rosenberg" Skin Graft Knives, Manual Dermatomes Unique Concept in Skin Grafting

New Generation of Ergonomic Manual Dermatomes

The "Rosenberg" Skin Graft Knives System is comprised of three dermatomes a Macro (standard), Mini & Micro design, suitable for all requirements and body areas. The design of these dermatomes incorporates an ergonomic approach that provide the surgeon with improved safety, control and ease of use. Its custom design facilitates tangential excisions & tissue harvesting in all body areas where traditional skin graft knives are used and is particularly suitable for the abdomen, chest, shoulders and back.

Unique Micrometer Adjustment

This unique depth regulation system incorporates the following advantages:

- Single micrometer adjustment.
- Guard always parallel with blade edge.
- Free access in restricted areas.

Practically maintenance-free

The "Rosenberg" Dermatome Features:

- Smoothly contoured one piece guard that protects the blade corners. This feature enables use of the entire length of the blade including the blade's corner that allows accurate tangential excision of irregular and difficult areas.
- The uniquely designed and positioned ergonomic handles maintain a clear operative field and provide for accurate, effortless control.
- The distinctive handle-blade design enables use on flat surfaces (abdomen, back), corners and other difficult areas.
- Simple, easy and safe loading and removal of standard blade.
- Positive blade-lock to prevent accidental opening during use.
- Superfine electro-polished guard that eliminates the need for mineral oil.



The Macro Dermatome

The standard (Macro) Skin Graft Knife may be used on all areas where traditional patterns of skin graft knives are employed (i.e. thighs), and is particularly suitable for the abdomen, chest, shoulders and back. The "Rosenberg" Macro Dermatome uses standard disposable skin graft blades.

The Mini Dermatome

The Mini Skin Graft Knife may be used on all areas where traditional patterns of skin graft knives are employed (i.e. thighs), and in flat areas where other dermatomes cannot be used. It may be used for a precise and controlled tangential excision (debridement and graft harvesting) of any shape and depth in any body area (i.e. hands, face & feet). Beside debridement and graft harvesting these features make it an ideal tool for a precise and controlled depithelialization surgery (i.e. mammoplasty) and preparation of dermal grafts.

The Micro Dermatome

This unique size (the smallest in the world) and slim shape of this Skin Graft Knife has been designed for very fine tangential excision and de-epidermization from areas such as fingers, nipples and face. It is also useful for small cutaneous lesion removal and graft harvesting from buccal mucosa, retro-auricular region, cartilage and perichondrium. It may be used for shaping and splitting cartilage grafts in reconstructive procedures.



Easy access during surgical intervention



Close-up of the Micro Dermatome



The Macro Dermatome



The Mini Dermatome




The Micro Dermatome



The "Rosenberg" Macro, Mini & Micro

Skin Graft Knives

New Generation of Ergonomic Manual Dermatomes

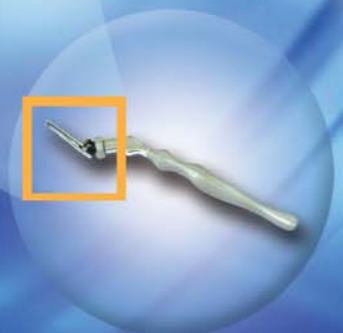

wallish adv.

Manufacturer and distributor:

• 4Med Ltd. Nahal-Oz, 85145 Israel • Tel: 972-8-6805223 • Fax: 972-8-6805231 • E-mail: info@4med.co.il
• www.4med.co.il

• 4Med US sales office: 63 Pond View Lane • Chappaqua NY 10514 USA • Phone: (914)238-8710
• Fax: Toll-free in USA and Canada (888)475-5089 • E-mail: sales@4med-surgical.com
• www.4med-surgical.com